2.2.7 Quasar Wireless Bluetooth Interface

Straight from your measuring sensor to your PC or Android device with no cables

- Quasar wireless interface connects to any Ophir sensor and broadcasts to your PC or your Android device running StarViewer
- Wireless range of 10-30 meters depending on surroundings
- Operates from rechargeable battery with typically >40 hours lifetime
- Powerful USB interface with StarLab PC application software included or StarViewer Android application
- Converts your PC or your android Device into a complete laser power/energy meter
- Log power and energy, average, statistics, histograms and more (only in PC)
- Monitor up to 7 Quasars simultaneously on one PC (only in PC)

Quasar Bluetooth Wireless Sensor to PC Interface

in e





Quasar module connects to any Ophir sensor, thermal, pyroelectric or photodiode

Any PC, laptop or Android device connects to Quasar module via Bluetooth adapter and operates as a power/energy meter/data logger

Specification


Sensor Compatibility	All Ophir standard sensors, Thermal ^(a) , Photodiode ^(b) and Pyroelectric							
Number of Sensors on One PC	Up to 7 Quasars can operate simultaneously and be displayed at the same time on one PC							
Operating Range	10-30 meters depending on surroundings when used with built in laptop Bluetooth or Ophir recommended adapter							
Power	Powered by rechargeable NiMH battery. Battery life typical 40 hours, 20 hours for pyro sensors. Automatically goes into sleep mode when not connected to PC. Low batt indication. Charges from 12VDC either polarity. The charger can be ordered from your local distributor.							
LED Indicator	LED indicator indicates whether connected, in standby or off							
Bluetooth Standard	Bluetooth class 1. Connection to PC is transparent to user. Will work with built in laptop Bluetooth and most add on USB to Bluetooth adapters.							
Data Transfer Rate for Pyro Sensors	500Hz							
Dimensions	94mm L x 96mm W x 36mm H not including antenna							
Connections	15 pin D type sensor connector standard Ophir 12V charger input							
Compliance	CE, UKCA, China RoHS							
Notes:	(a) When operating with BeamTrack sensors, measures Power & Energy only. (b) Not including BC20, PD300-CIE and PD300RM sensors.							

Ordering Information

Item	Description	Ophir P/N
Quasar Bluetooth Interface	Module to operate one Ophir sensor from your PC via Bluetooth wireless interface. Comes with software. Max repetition rate for every pulse 500Hz. Powered from built in rechargeable battery. Comes with power supply. Bluetooth adapter required when not available on PC. See next line	7Z01300
Battery Pack for Quasar	Replacement battery pack for Quasar	7E14007A
N Polarity power supply/charger	Power Supply/Charger AC/DC 12V 2A N-2.1x5.5 (1 unit supplied with Quasar)	7E05029



2.2 PC Interfaces 2.2.1 PC Connectivity Options for Power/Energy Measurement

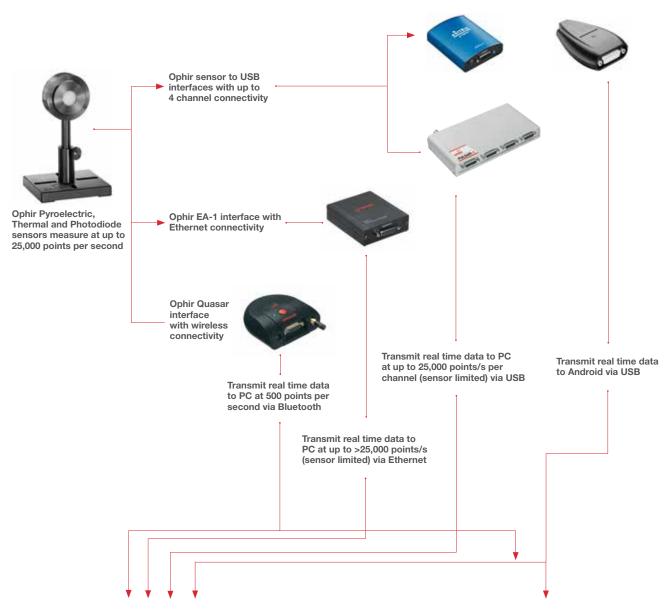

StarLab Software (data transmitted via USB, Ethernet or Bluetooth)



StarLab Software

StarViewer Application (data transmitted via Bluetooth and USB)





StarViewer Android Application

Power Meters

For latest updates, please visit our website: www.ophiropt.com

2.2.8 Summary of Computer Options for Ophir Meters and Interfaces

Communications

With Ophir RS232, GPIB, Bluetooth, USB and Ethernet communication options you can transfer data from the sensor to the computer in real time or offline. You can also control your Ophir power meter from the computer.

- USB on Nova II, Vega, StarBright, Centauri (optional on StarLite) power meters and Juno, Juno+, Pulsar PC interfaces
- Bluetooth wireless on Quasar interface
- RS232 on LaserStar, Nova II, Vega, StarBright, Centauri and Juno-RS optional on Nova
- GPIB optional on LaserStar
- Ethernet on EA-1 interface

Ophir Power Meter and Interface Specifications

• p • .			acc opec	inoutione			_				
Model	Centauri	StarBright	Nova II / Vega	StarLite	LaserStar	Nova	Juno / Juno+	Juno-RS	Pulsar-1, 2 or 4	EA-1	Quasar Bluetooth
Communication method	USB / RS232	USB / RS232	USB / RS232	USB (c)	RS232 / GPIB	RS232	USB	RS232	USB	Ethernet	Bluetooth
Power Measuren	nent								-		
Power log period	1s to 1000hr.	1s to 1000hr.	12s to 600hr.	N.A	12s to 600hr.	5s to 24hr.	1s to Unlimited	1s to Unlimited	1s to Unlimited	1s to Unlimited	1s to Unlimited
Max points stored onboard	Unlimited	Unlimited	Nova II 5400 Vega 27000	N.A	5400	300	N.A	N.A	N.A	N.A	N.A
Max points direct on PC	Unlimited	Unlimited	Unlimited	N.A	Unlimited	Unlimited	Unlimited	Unlimited	Unlimited	Unlimited	Unlimited
Analog output	1V, 2V, 5V, 10V F.S.	1V, 2V, 5V, 10V F.S.	1V, 2V, 5V, 10V F.S.	1V F.S.	1V F.S.	1V F.S.	N.A / 1V, 2V, 5V, 10V F.S.	1V, 2V, 5V, 10V	N.A	N.A	N.A
Energy Measurer	nent										
Max logging rate	25,000Hz USB 30Hz RS232	5000Hz USB 30Hz RS232	>2000Hz USB ^(a) >30Hz RS232	20Hz ^(c)	>30Hz RS232 >1500Hz GPIB ^(a)	>10Hz	10,000Hz ^(a)	500Hz ^(a)	25,000Hz ^(a)	>25,000Hz ^(a)	500Hz
Max onboard data logging rate	25,000Hz	5000Hz	4000Hz ^(a)	N.A	>1500Hz ^(a)	>10Hz	N.A	N.A	N.A	N.A	N.A
Max points stored USB/onboard	Unlimited	Unlimited	Nova II 59,400 Vega 250,000	N.A	59,400	1000	N.A	N.A	N.A	N.A	N.A
Trigger input and output	Trigger input to synchronize measurement of pulses	N.A	N.A	N.A	N.A	N.A	N.A	N.A	BNC trigger input to enable measurement of missing pulses. Can also be configured to give trigger output	N.A	N.A
Timing - time stamp for each pulse	resolution 1µs	resolution 1µs	N.A	N.A	N.A	N.A	resolution 1µs	resolution 1µs	resolution 1µs	resolution 1µs	resolution 10ms
General											
Com Object	yes	yes	yes	yes ^(c)	no	no	yes	no	yes	yes	no
LabVIEW VIs Maximum baud rate	yes 115200	yes 115200	yes 38400	yes ^(c) N.A	yes 38400	yes 19200 ^(b)	yes N.A.	no 115200	yes N.A.	N.A.	no N.A.
PC file format					Text files, sprea	debeet eemee					
TTL Out	ves	N.A	N.A	N.A	N.A	N.A	N.A	N.A	N.A	N.A	N.A
Number of sensors supported	2 / 1 sensors per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC		One sensor per unit	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit	4 / 2 / 1 sensors per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Car combine several units with softwar for display of up to 7 Quasars on one PC
Compatible sensors				Supports me	ost Ophir pyroele	ectric, thermal a	and photodioc	le sensors			
Power supply	Powered from internal rechargeable battery power supply	Powered from internal rechargeable battery power supply	Powered from internal rechargeable battery power supply	Powered from internal rechargeable battery power supply	Powered from internal rechargeable battery power supply	Powered from internal rechargeable battery power supply	Powered from USB	12V wall cube plugs into jack on rear	12V wall cube plugs into jack on rear	12V wall cube plugs into jack or PoE	Powered from internal rechargeable battery powe supply
Dimensions	47 x 200 x 130mm	212 x 114 x 40mm	208 x 110 x 43mm / 210 x 109 x 36mm	211 x 114 x 40mm	194 x 228 x 57mm	205 x 95 x 39mm	77 x 55 x 23mm / 105 x 80 x 29mm	114 x 80 x 29mm	103 x 190 x 33mm	93 x 73 x 29mm	94 x 96 x 36mm
Notes:	 (b) For pyroelectric (c) StarLite must b 	c sensors, maxim e USB enabled i	logging every single num guaranteed ba n order to work with enabled (P/N 7Z015	ud rate is 9600. n StarLab. If your							Activation Code

2.3 Software Solutions 2.3.1 StarLab

StarLab turns your PC into a laser power/energy multi-channel station

Extensive Graphic Display of Data

- Line Plot, Histogram, Bar chart, Simulated Analog Needle
- Multiple data sets on one graph or separate graphs on the same screen

Advanced Measurement Processing

- Power/Energy Density, Scale Factor, Normalize against a reference
- Multi-channel comparisons
- User defined mathematical equations: channels A/B, (A-B)/C etc.
- Position & size measurement with BeamTrack sensors

Flexible Display Options with StarLab

Data Logging for Future Review


- Can be displayed graphically or saved in text format
- Easily exported to an Excel spreadsheet

Fully supports IPM, Ariel, Centauri, StarBright, StarLite, Vega, Nova II, Pulsar, Juno, Juno+, Juno-RS, Quasar and EA-1 devices with all standard Ophir sensors


You may choose to display them separately

Choose which channels to display Maximize one of the sources 1.967uW 75.40mW Ø 84.10 477.6ell 1.970 3A-P-V1 (s/n: 619524) 50(150)A-PPS (3/h: 643979) 1.217vW \$46.049 Vega (s/n: 570904) Juno (s/n: 345003) 2.668uJ Ø 32.79uW PD300-JRG (s/n: 105579) PD300 fs/m: 2 PE10-C (s/n: 333010) Open sensors in new window Setup screen Choose line graph 2 ----曰 130.7mW 110.7mW 110.7mW : 1 = 0 or histogram 0.12 0.18 One of the above screens is maximized 117.4mW Min Max 188.2mW 0.000mW or needle display


Multiple Sensors displayed together

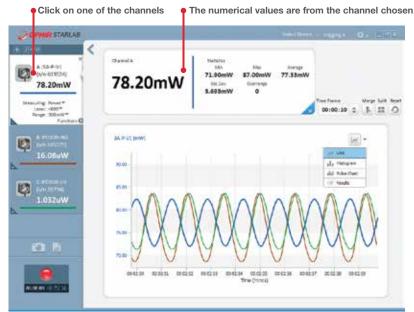

Here multi line graph display has been chosen

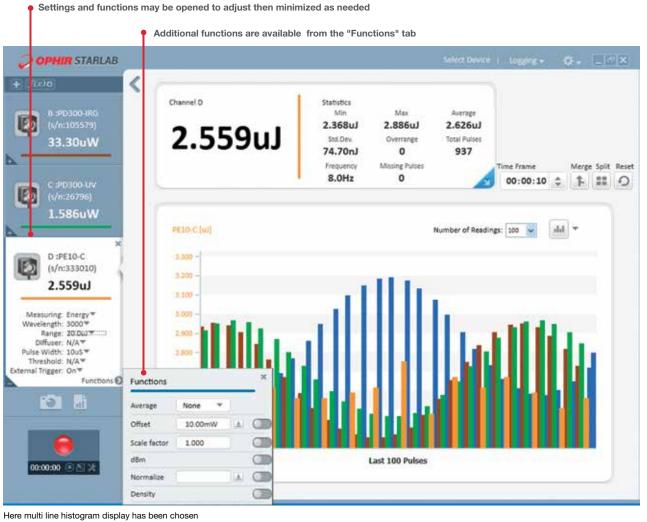



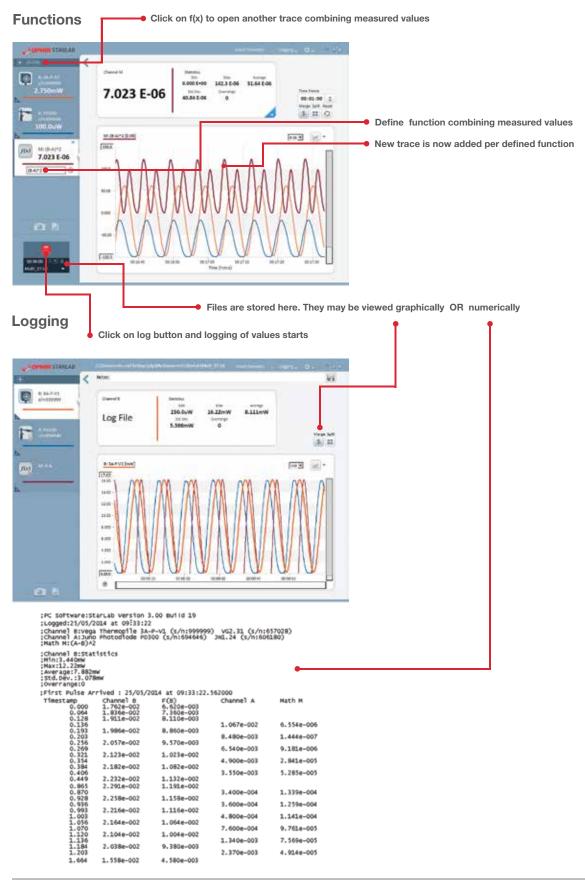

BeamTrack Power/Position/Size Screens

Open Measuring type tab and choose Track

OPHIR STARLAB < Channel A Statistics Min A: 50(150)A-PPS Max Average ō s/n:643979 5.680W 5.700W 5.684W 5.690W 5.690W Overrange Std.Dev • Power 5.072mW 0 Measuring: Track w/ Pov Laser: 10.6 * Range: 50.0W Power Energy Func ns O Track w/ Pov A: 50(150)A-PPS[mm] 0 * 10 -1.07mm X: Position +0.44mm 2.92mm Size: -10 10 191 M Size 00:00:00 (C) X -10

Power / Position / Size screen

