

2.2.4 Juno-RS Interface

Virtual meter (PC-Interface) allowing connectivity of Ophir's smart (DB-15) power and energy sensors to automation systems and computers using RS-232 communication

- Plug and play with all standard Ophir smart sensors
- Analog output in autonomous mode: Outputs voltage relative to measurement while connected standalone to a sensor
- Robust connectors, suited for semi-industrial environment
- Advanced logging and data processing with included StarLab application
- Pulsed Power measurements with Thermopile detectors
- Low Frequency Power - power measurement from pulse cycle energy (for VCSEL)




Smart Sensor to Juno-RS to PC

Ophir's Juno-RS, virtual meter provides easy integration into automation systems (e.g. PLC) using RS-232 in their communication infrastructure. It also allows turning the customer's PC or laptop into a full-fledged Ophir laser power/energy meter. The sensor is connected to the Juno-RS which sends the measurement results to the PLC or the destination computer.

The Juno-RS has an analog output which provides an output voltage proportional to the power/energy measured. When using the Juno-RS with a PC/Laptop, just install the StarLab software, plug the sensor into the Juno-RS module and connect the Juno-RS with a standard RS-232 cable to a PC serial port or USB to RS-232 adapter. You can connect several Juno-RS modules to the PC at the same time.



RS-232 Terminal

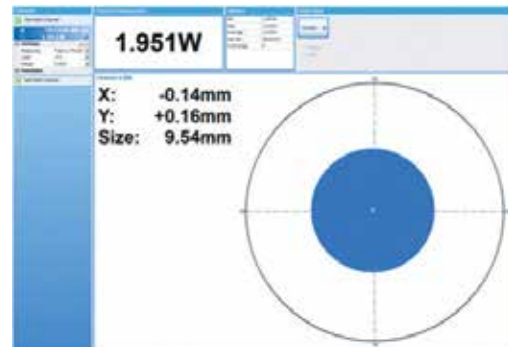

Juno-RS with BeamTrack sensor and StarLab showing beam power, position and size

Specifications


Power Measurement	
Power log period	Unlimited
Energy Measurement	
Max logging rate	500Hz @ baud rate 115200 (a)
Trigger input and output	N.A.
Timing	Supports time stamp for each pulse - resolution 1µs
General	
Number of sensors supported	One sensor per unit. Can combine several units with StarLab software for display of up to 8 sensors on one PC
Compatible sensors	Supports all standard Ophir Pyroelectric (PE-C series), Thermal, BeamTrack and Photodiode sensors. Works with our PD300RM sensors
Power supply	Industrial +12V wall cube power supply, plugs into jack on rear. The power supply is supplied with the device or can be ordered from your local distributor
Power consumption	100mA max at 12V
Power input voltage	11V min to 25V max (if using alternative to provided power supply)
Outputs	RS-232; Analog output with user selectable full-scale of 1, 2, 5 or 10V
Dimensions	114mm L x 80mm W x 29mm H
Weight	0.174 kg
Compliance	CE, UKCA, China RoHS
Note: (a) This is the data logging rate for every single point in turbo mode. Above that rate, the instrument will sample points but not log every single point	

Ordering Information

Item	Description	Ophir P/N
Juno-RS	Module to operate one Ophir sensor from your PC RS-232 port. Comes with software	7Z01254
Juno-RS RS-232 Cable	D9 Male/Female 1.8-2 meter (1 unit supplied with Juno-RS)	7E11216
Juno-RS Power Supply	Power Supply AC/DC 12V 2A 2.5x5.5x13.8 S (1 unit supplied with Juno-RS)	7E05093
Juno-RS Analog Output Cable	Coax BNC-M to SMA-M RG-174 Cable 2 meter	7E01541

2.2 PC Interfaces

2.2.1 PC Connectivity Options for Power/Energy Measurement

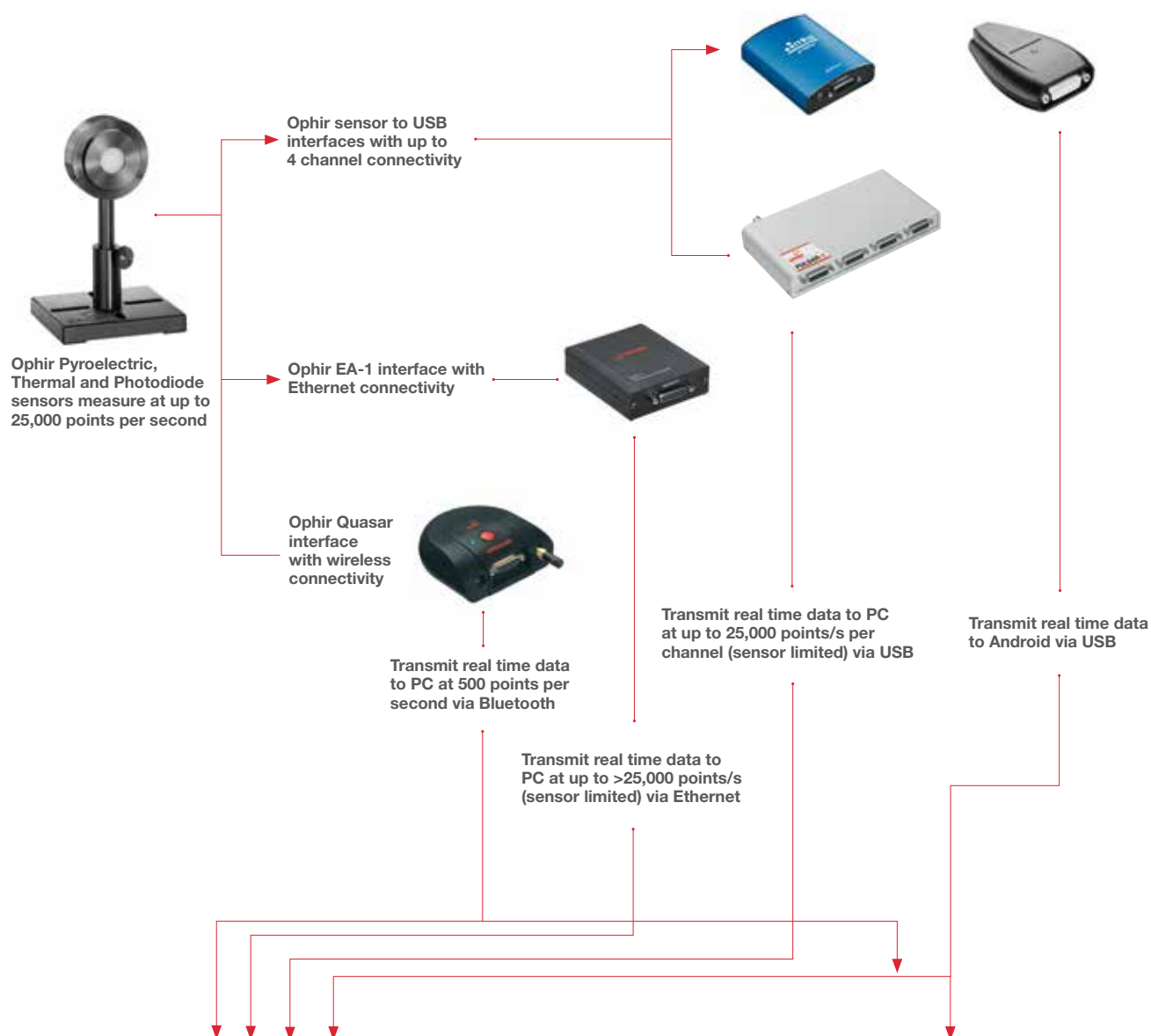

StarLab Software (data transmitted via USB, Ethernet or Bluetooth)

StarViewer Application (data transmitted via Bluetooth and USB)



StarLab Software



StarViewer Android Application

2.2.8 Summary of Computer Options for Ophir Meters and Interfaces

Communications

With Ophir RS232, Bluetooth, USB and Ethernet communication options you can transfer data from the sensor to the computer in real time or offline. You can also control your Ophir power meter from the computer.

- USB on Nova II, Vega, StarBright, Centauri (optional on StarLite) power meters and Juno, Juno+, Pulsar PC interfaces
- Bluetooth wireless on Quasar interface
- RS232 on Nova II, Vega, StarBright, Centauri and Juno-RS optional on Nova
- Ethernet on EA-1 interface and Centauri power meter

Ophir Power Meter and Interface Specifications

Model	Centauri	StarBright	Nova II / Vega	StarLite	Nova	Juno / Juno+	Juno-RS	Pulsar-1, 2 or 4	EA-1	Quasar Bluetooth
Communication method	USB / RS232 / Ethernet	USB / RS232	USB / RS232	USB ^(c)	RS232	USB	RS232	USB	Ethernet	Bluetooth
Power Measurement										
Power log period	1s to 1000hr.	1s to 1000hr.	12s to 600hr.	N.A	5s to 24hr.	1s to Unlimited	1s to Unlimited	1s to Unlimited	1s to Unlimited	1s to Unlimited
Max points stored onboard	Unlimited	Unlimited	Nova II 5400 Vega 27000	N.A	300	N.A	N.A	N.A	N.A	N.A
Max points direct on PC	Unlimited	Unlimited	Unlimited	N.A	Unlimited	Unlimited	Unlimited	Unlimited	Unlimited	Unlimited
Analog output	1V, 2V, 5V, 10V F.S.	1V, 2V, 5V, 10V F.S.	1V, 2V, 5V, 10V F.S.	1V F.S.	1V F.S.	N.A / 1V, 2V, 5V, 10V F.S.	1V, 2V, 5V, 10V	N.A	N.A	N.A
Energy Measurement										
Max logging rate	25,000Hz USB 30Hz RS232	5000Hz USB 30Hz RS232	>2000Hz USB ^(a) >30Hz RS232	20Hz ^(c)	>10Hz	10,000Hz ^(a)	500Hz ^(a)	25,000Hz ^(a)	>25,000Hz ^(a)	500Hz
Max onboard data logging rate	25,000Hz	5000Hz	4000Hz ^(a)	N.A	>10Hz	N.A	N.A	N.A	N.A	N.A
Max points stored USB/onboard	Unlimited	Unlimited	Nova II 59,400 Vega 250,000	N.A	1000	N.A	N.A	N.A	N.A	N.A
Trigger input and output	Trigger input to synchronize measurement of pulses	N.A	N.A	N.A	N.A	N.A	N.A	BNC trigger input to enable measurement of missing pulses. Can also be configured to give trigger output	N.A	N.A
Timing - time stamp for each pulse	resolution 1µs	resolution 1µs	N.A	N.A	N.A	resolution 1µs	resolution 1µs	resolution 1µs	resolution 1µs	resolution 10ms
General										
Com Object	yes	yes	yes	yes ^(c)	no	yes	no	yes	yes	no
LabVIEW VIs	yes	yes	yes	yes ^(c)	yes	yes	no	yes	no	no
Maximum baud rate	115200	115200	38400	N.A	19200 ^(b)	N.A.	115200	N.A.	N.A.	N.A.
PC file format										
TTL Out	yes	N.A	N.A	N.A	N.A	N.A	N.A	N.A	N.A	N.A
Number of sensors supported	2 / 1 sensors per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC	4 / 2 / 1 sensors per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 8 sensors on one PC	One sensor per unit. Can combine several units with software for display of up to 7 Quasars on one PC
Compatible sensors	Supports most Ophir pyroelectric, thermal and photodiode sensors									
Power supply	Powered from internal rechargeable battery power supply	Powered from internal rechargeable battery power supply	Powered from internal rechargeable battery power supply	Powered from internal rechargeable battery power supply	Powered from internal rechargeable battery power supply	Powered from USB	12V wall cube plugs into jack on rear	12V wall cube plugs into jack on rear	12V wall cube plugs into jack or PoE	Powered from internal rechargeable battery power supply
Dimensions	47 x 200 x 130mm	212 x 114 x 40mm	208 x 110 x 43mm / 210 x 109 x 36mm	211 x 114 x 40mm	205 x 95 x 39mm	77 x 55 x 23mm / 105 x 80 x 29mm	114 x 80 x 29mm	103 x 190 x 33mm	93 x 73 x 29mm	94 x 96 x 36mm

Notes: (a) The above refers to the rate for logging every single point in turbo mode. Above that rate, the instrument will sample points but not log every single point.
(b) For pyroelectric sensors, maximum guaranteed baud rate is 9600.
(c) StarLite must be USB enabled in order to work with StarLab. If your StarLite has not been USB enabled, please contact your Ophir distributor in order to obtain a USB Activation Code. Now available: StarLite with USB enabled (P/N 7Z01569)

2.3 Software Solutions

2.3.1 StarLab

StarLab turns your PC into a laser power/energy multi-channel station

Extensive Graphic Display of Data

- Line Plot, Histogram, Bar chart, Simulated Analog Needle
- Multiple data sets on one graph or separate graphs on the same screen

Advanced Measurement Processing

- Power/Energy Density, Scale Factor, Normalize against a reference
- Multi-channel comparisons
- User defined mathematical equations: channels A/B, (A-B)/C etc.
- Position & size measurement with BeamTrack sensors

Data Logging for Future Review

- Can be displayed graphically or saved in text format
- Easily exported to an Excel spreadsheet

Fully supports IPM, Ariel, Centauri, StarBright, StarLite, Vega, Nova II, Pulsar, Juno, Juno+, Juno-RS, Quasar and EA-1 devices with all standard Ophir sensors

Flexible Display Options with StarLab

Choose which channels to display



Setup screen




One of the above screens is maximized

You may choose to display them separately


Maximize one of the sources



Choose line graph

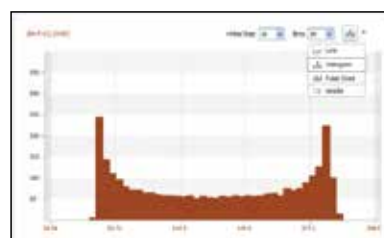

or histogram

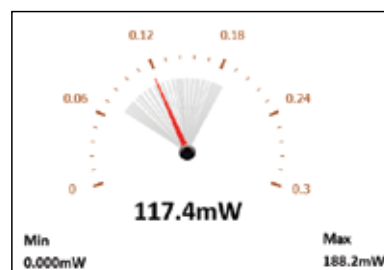

or needle display

Multiple Sensors displayed together

- Click on one of the channels
- The numerical values are from the channel chosen



Here multi line graph display has been chosen

- Settings and functions may be opened to adjust then minimized as needed
- Additional functions are available from the "Functions" tab



Here multi line histogram display has been chosen

Functions and Logging

Functions

Click on f(x) to open another trace combining measured values



Define function combining measured values

New trace is now added per defined function

Logging

Click on log button and logging of values starts

Files are stored here. They may be viewed graphically OR numerically



```


:PC Software:StarLab version 3.00 Build 19
:Logged:25/05/2014 at 09:33:22
:Channel B:Vega Thermopile 3A-P-V1 (s/n:999999) VG2.31 (s/n:657028)
:Channel A:Juno Photodiode P0300 (s/n:694646) JN1.24 (s/n:606180)
:Math M:(A-B)*2
:Channel B:Statistics
:Min:3.440mW
:Max:12.22mW
:Average:7.882mW
:Std.Dev.:3.078mW
:Overrange:0
:First Pulse Arrived : 25/05/2014 at 09:33:22.562000

```


Timestamp	Channel B	F(B)	Channel A	Math M
0.000	1.762e-002	6.620e-003		
0.064	1.836e-002	7.350e-003		
0.128	1.911e-002	8.110e-003		
0.136			1.067e-002	6.554e-006
0.193	1.986e-002	8.860e-003		
0.203			8.480e-003	1.444e-007
0.256	2.057e-002	9.570e-003		
0.269			6.540e-003	9.181e-006
0.321	2.123e-002	1.023e-002		
0.354			4.900e-003	2.841e-005
0.384	2.182e-002	1.082e-002		
0.406			3.550e-003	5.285e-005
0.449	2.232e-002	1.132e-002		
0.465				
0.470	2.291e-002	1.191e-002		
0.870			3.400e-004	1.339e-004
0.928	2.258e-002	1.158e-002		
0.936			3.600e-004	1.259e-004
0.993	2.216e-002	1.116e-002		
1.003			4.800e-004	1.141e-004
1.056	2.164e-002	1.064e-002		
1.070			7.600e-004	9.761e-005
1.120	2.104e-002	1.004e-002		
1.136			1.340e-003	7.569e-005
1.184	2.038e-002	9.380e-003		
1.203			2.370e-003	4.914e-005
1.664	1.558e-002	4.580e-003		

BeamTrack Power/Position/Size Screens

- Open Measuring type tab and choose Track

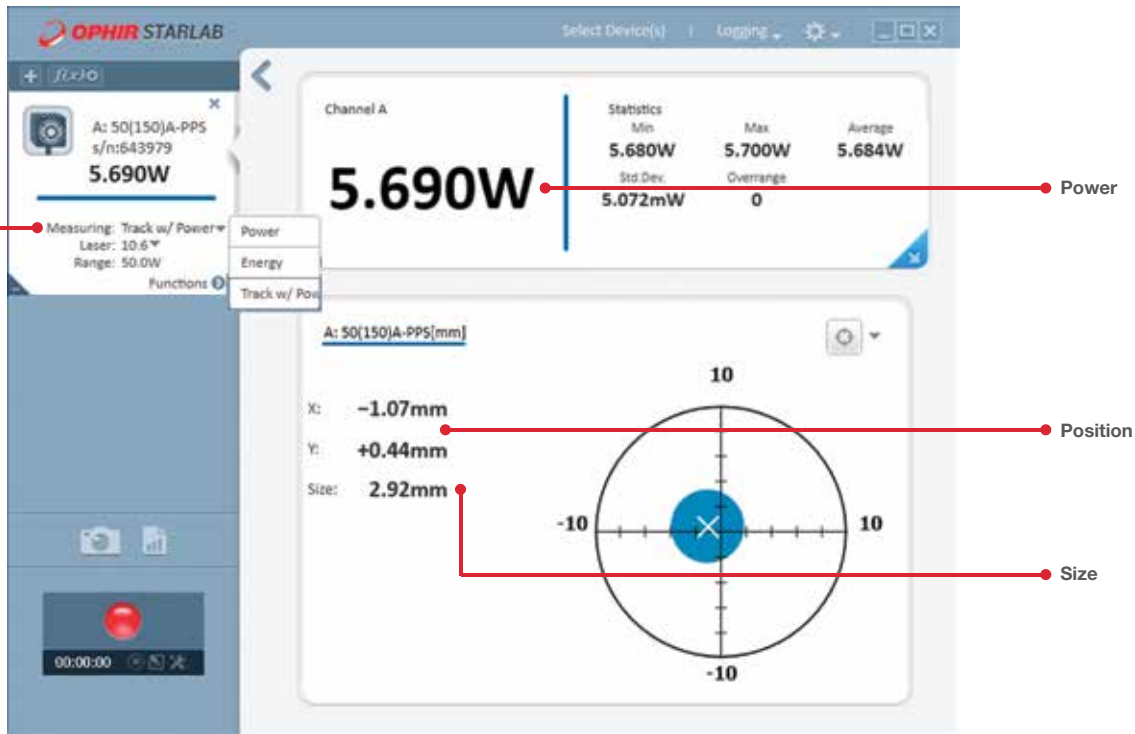

Power / Position / Size screen

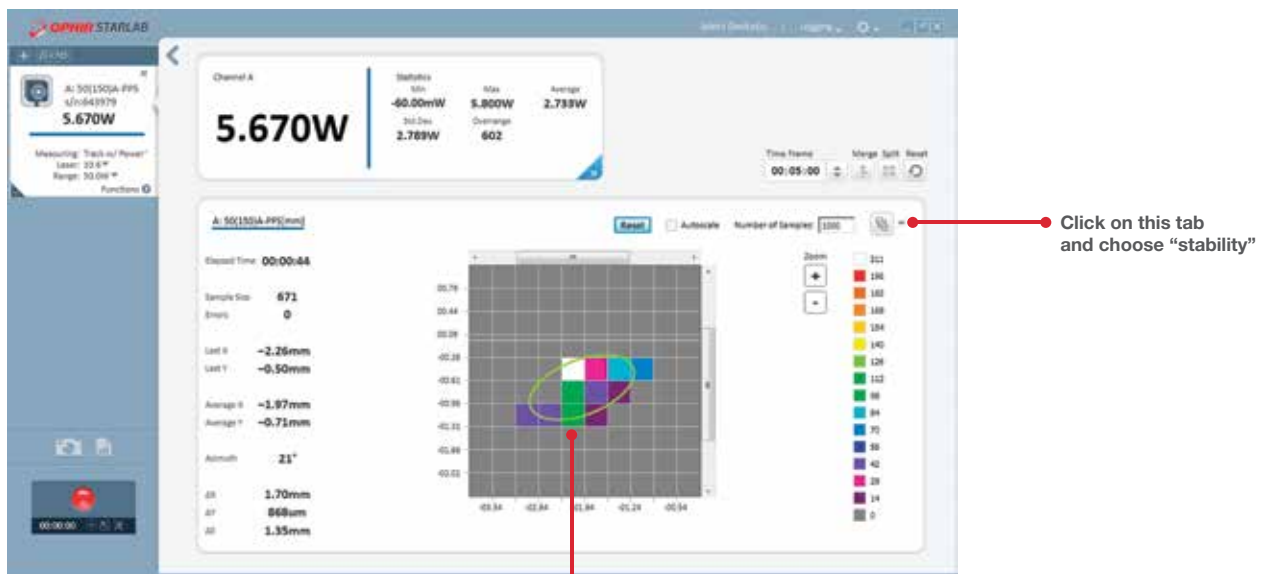

Position stability screen

- Displays beam center wander weighted for dwell time at each position