

1.1.1.3 Special Photodiode Sensors

50pW to 50mW

Features

- PD300-BB for broadband light sources - radiometry (PD300-BB-50mW option up to 50mW)
- BC20 for measuring scanned beams such as bar code light sources

PD300-BB / PD300-BB-50mW



BC20



Model	PD300-BB	PD300-BB-50mW		BC20
Use	Radiometry-broad spectrum	Same as PD300-BB with removable attenuator for use to 50mW		Scanned beams e.g. bar code with continuous wavelength curve
Detector Type	Silicon with special filter	Silicon with special filter		Silicon with peak and hold circuit
Aperture	10x10mm	10x10mm		10x10mm
Spectral Range nm	430 - 1000 (see graph)	430 - 1000 (see graph)		400 - 1100 (see graph) ^(b)
Calibration Uncertainty ±%	1.1 430-1000nm ^(c)	1.1 430-1000nm ^(c)		NA
Filter Mode		Filter out	Filter in	
Power Range	50pW to 4mW	50pW to 4mW	1nW to 50mW	0.1mW to 20mW
Power Scales	4mW to 8nW and dBm	4mW to 8nW and dBm	50mW to 80nW and dBm	20mW to 2mW
Resolution	1pW	1pW	10pW	1μW
Accuracy	Maximum deviation from flat spectrum (see graph) ±10%	Maximum deviation from flat spectrum (see graph) ±10%	430-910nm, ±12% flatness	±3% for >10% of full scale. Deviation from calibration -3% at 30,000 inch/s scan rate on sensor
Damage Threshold W/cm ²	10	10	100	50
Max Pulse Energy μJ	1	1	10	NA
Noise Level pW	2	2	30	5μW
Response Time with Meter s	0.2	0.2	0.2	Two modes of operation: Hold: holds highest reading for 5s then updates. No Hold: updates reading 3 times per second
Beam Position Dependence	±2% for broadband light sources	±2% for broadband light sources	±3% for broadband light sources	±2%
Background Subtraction	NA	NA	NA	Background is automatically subtracted from both scanned and static beams
Fiber Adapters Available (see page 35)	NA	ST, FC, SMA, SC		NA
Compatible Meter /Interface	All Meters & Interfaces	All Meters & Interfaces		StarBright, Vega, Nova II, Juno, Juno+, LaserStar and Nova
Compliance	CE, UKCA, China RoHS	CE, UKCA, China RoHS		CE, UKCA, China RoHS
Version				V1
Part Number	7Z02405	7Z02440		7Z02481 ^(a)

Notes:


(a) Swivel stand for BC20 sensor P/N 1Z09004

(b) The user can select up to 5 wavelengths from the spectral range. When used with the Nova or LaserStar meters, the sensor will only have the discrete wavelengths 405nm, 633nm, 650nm, 675nm and 780nm


(c) For calibration uncertainty of wavelengths outside of this range see table on page 25

* For graphs see page 33-34


PD300-BB / PD300-BB-50mW with filter off

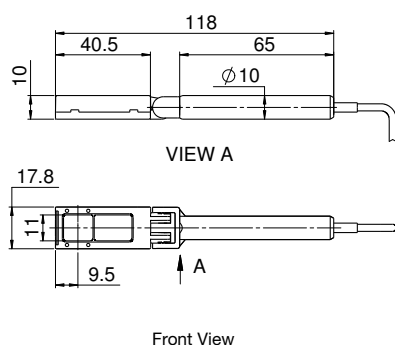

PD300-BB-50mW with filter installed

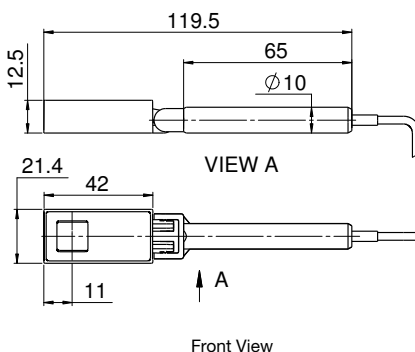


BC20

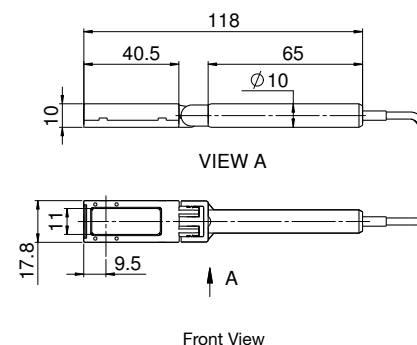


1.1.1.4 Graphs


Temperature Coefficient of Sensitivity

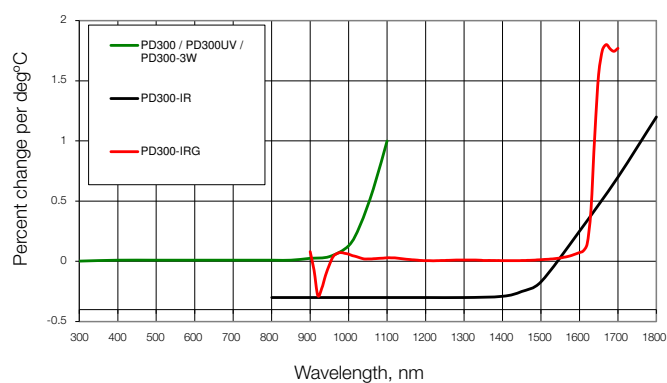

PD300 Angle Dependence

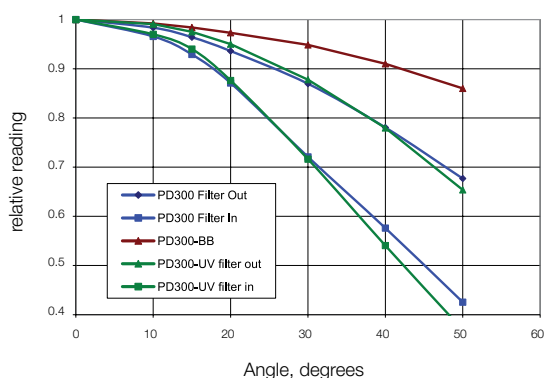


Dependence of Sensitivity on Numerical Aperture (PD300 - IRG)

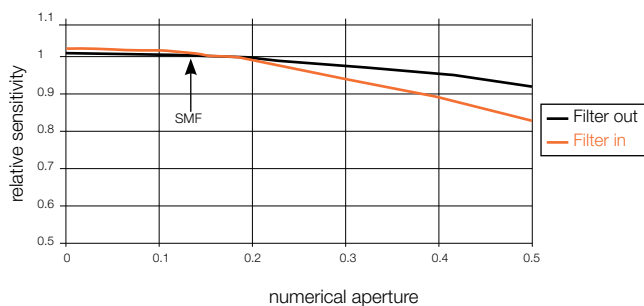


Note:

1. Graph assumes equal intensity into all angles up to maximum N.A.
2. Calibration is done with SMF, N.A. 0.13

Typical Sensitivity Curve of PD300-BB Sensors

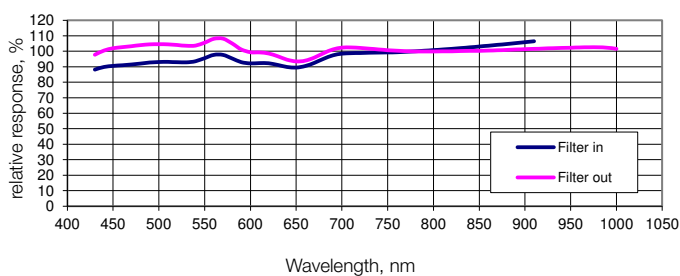


BC20 Relative Spectral Response

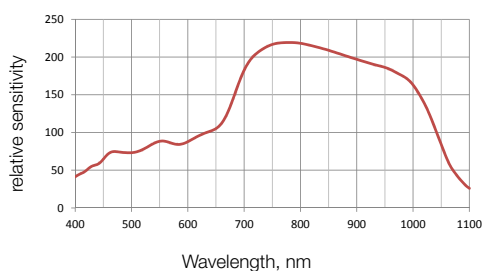


Approximate Spectral Response

Relative to 633nm or 1550nm


PD300 / PD300R

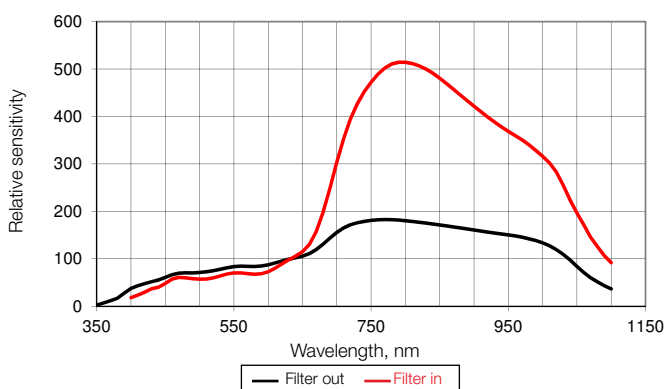


PD300-UV / PD300R-UV

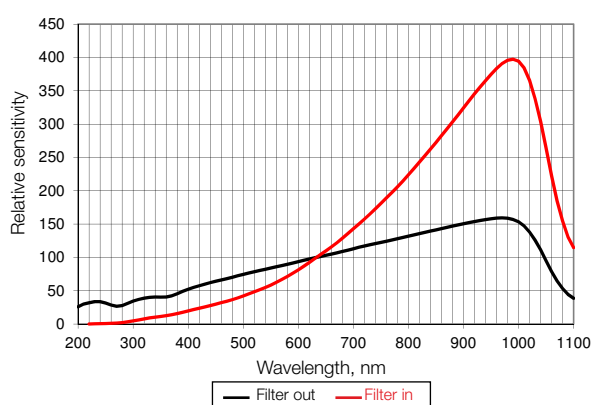


PD300-IR / PD300R-IR

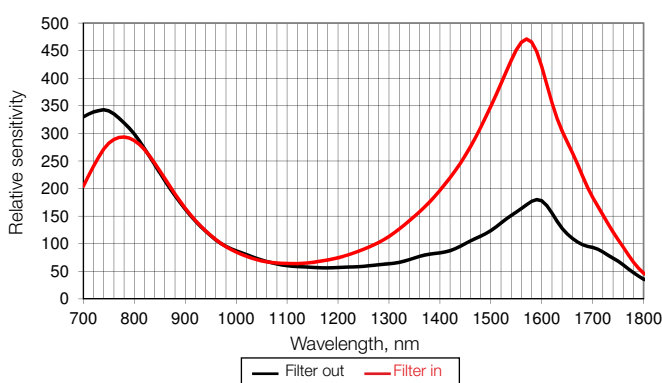


PD300-TP

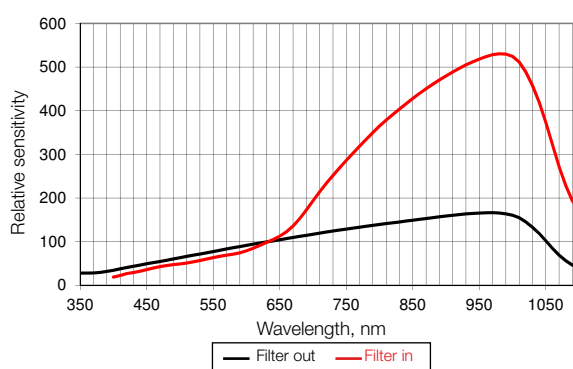

PD300-3W / PD300R-3W

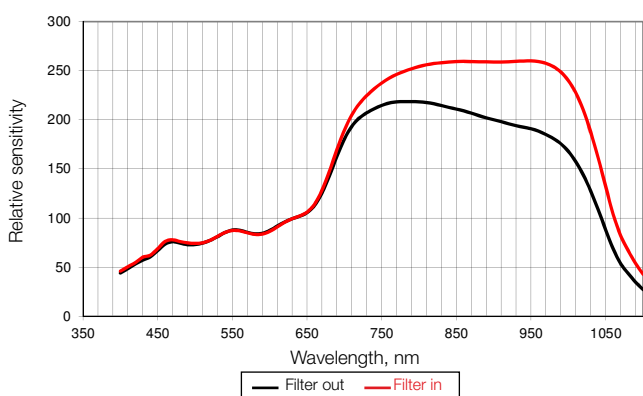

PD300-IRG

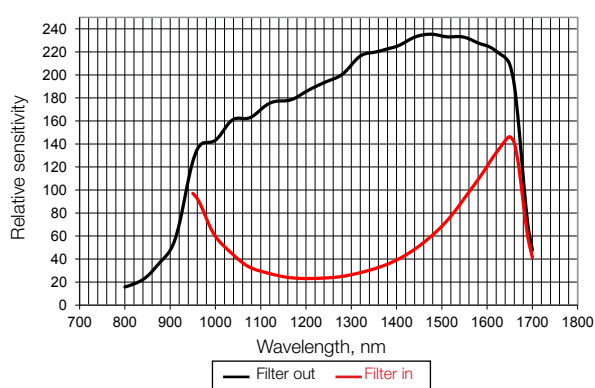

PD300-1W

